




CT AEC

Automatic Exposure Control in CT

CT AEC is used for constancy testing of CT AEC systems.

CT AEC includes:

- 11x PMMA ellipses of 25mm thickness and different height and width which can be re-configured to any order (ranging from the smallest 117 x 175.5mm to the largest 285.8 x 428.7mm)
- Each step contains a central hole to accommodate a choice of
 - Blank inserts for AEC checks (supplied)
 - Tubes to accommodate pencil ion chambers for dosimetry (supplied)
 - SedentexCT IQ image quality inserts (available separately)
- 2 x sturdy end plates with integral carry handles and adjustable height and length nylon supports.
- Spirit level and pigment-filled centre markings for easy and reproducible positioning.
- Wheeled protective storm case with telescopic handle.

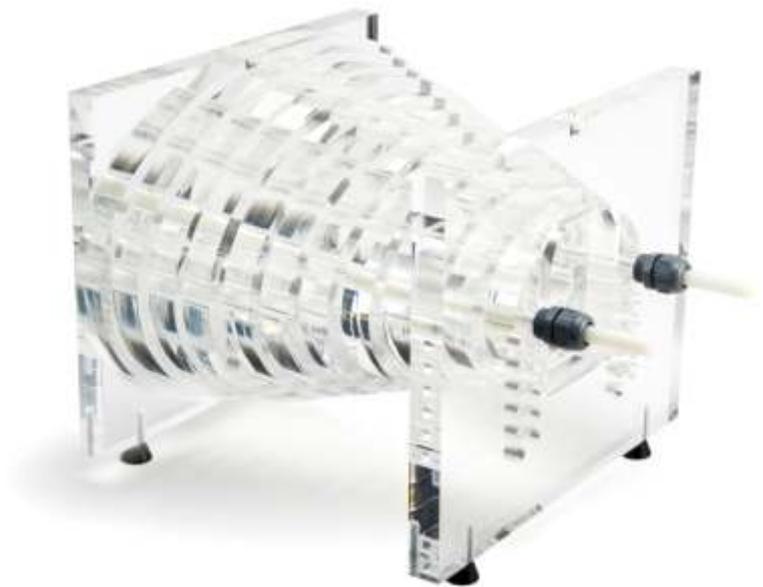

**FREE
Software
Included!**



CT AEC

CT

Product Detail



www.leadstestobjects.com